



MetaboCensoR Tutorial

2026-07-06

Tutorial Edition 1.3

Software Release 1.3.x



Contents

General Information.....	3
Running & Installation.....	4
Input Requirements	4
Overview of Interface	5
Implementation Details.....	6
Quick Start	6
Tab “Upload Data”: Control Panel	7
Tab “Upload Data”: Display Area	8
Tab “Blank Filters”: Control Panel	9
Tab “Blank Filters”: Display Area	10
Tab “MS Filters”: Control Panel.....	11
Tab “MS Filters”: Display Area.....	12
Tab “QC Filters”: Control Panel.....	13
Tab “QC Filters”: Display Area	14
Tab “Peak Filters”: Control Panel.....	15
Tab “Peak Filters”: Display Area.....	16
Tab “Final Summary”: Control Panel.....	17
Tab “Final Summary”: Display Area.....	18
Tab “About” Overview	19



General Information


This tutorial aims at guiding the user to process data LC-MS data through Shiny GUI interface of MetaboCensoR. MetaboCensoR is an input-versatile Shiny application and local R package for analyte-centric peak table filtering. The workflow integrates four complementary modules for blank filtering, redundant ion-species filtering, quality-control filtering, and peak-based filtering. MetaboCensoR also provides interactive threshold optimization, exportable annotation tables, and synchronized filtering of associated .mgf files.

An overview of MetaboCensoR is presented below. The workflow is organized into six functional tabs that collectively support peak table curation, parameter optimization, and export of filtered results. The output of the app is a filtered peak table and, optionally, an MGF spectra file. Only the peak table upload is required for full workflow execution and final output generation. All other tabs are optional.

License: GPL-3.0

Contact: plyushchenko.ivan@gmail.com

Project Repository: github.com/plyush1993/MetaboCensoR

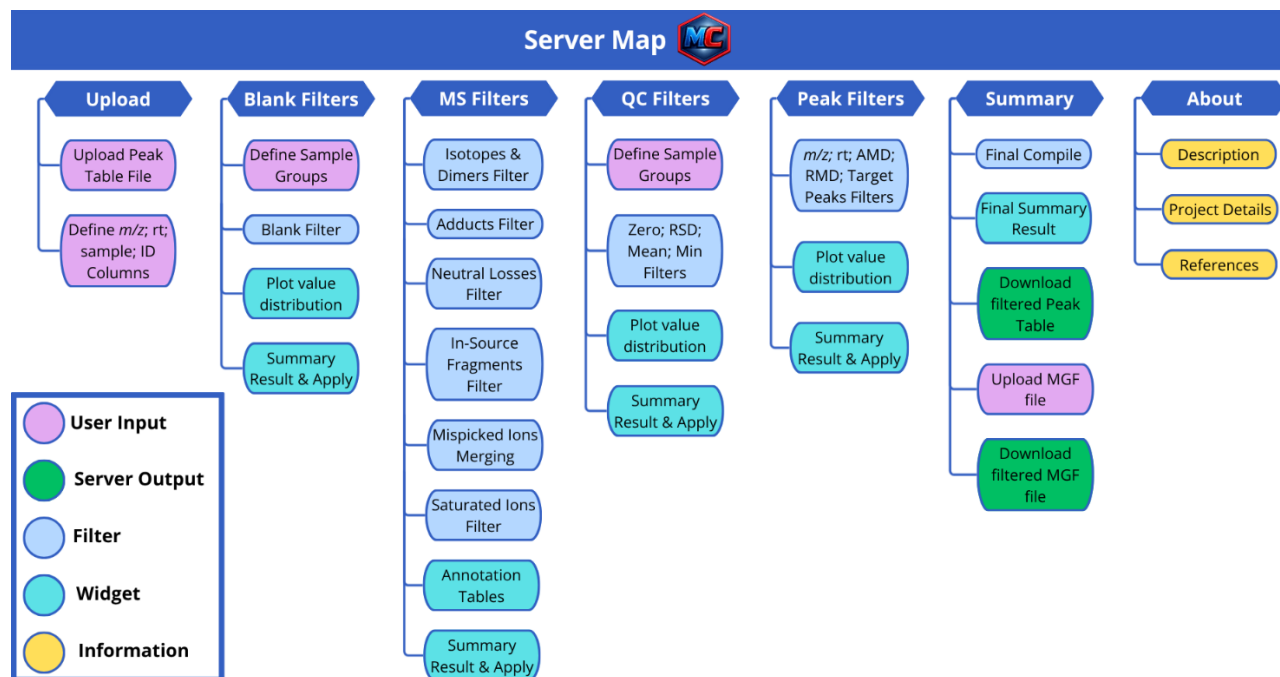



Running & Installation

Shiny Web Deployment is available by the link:

plyush1993.shinyapps.io/metabocensor/

Local installation of R package:

```
if (!require("BiocManager", quietly = TRUE)) {  
  install.packages("BiocManager")  
}  
if (!require("remotes", quietly = TRUE)) {  
  install.packages("remotes")  
}  
remotes::install_github("plyush1993/MetaboCensor", INSTALL_opts = "--no-multiarch")
```

or

```
if (!requireNamespace("pak", quietly = TRUE)) install.packages("pak")  
pak::pak("plyush1993/MetaboCensor")
```

Run Locally:

```
MetaboCensor::run_metabocensor()
```

Input Requirements

The application accepts peak tables in **.csv** format and MS1/MS2 mass spectral data in **.mgf** format, both of which can be directly exported from dedicated metabolomics processing software (xcms, MZmine, MS-DIAL). Make sure that your input peak table is aligned with data formats. Examples of acceptable formats are available here: [Input Examples](#).

Important note

Although the web version has been successfully tested on datasets with over 10,000 peaks, the Shiny server has limited memory capacity. Computationally intensive steps (like MS filters) or large MGF files (> 3 MB) may cause the app to stack or disconnect. For heavy processing, you should run the application locally.



Overview of Interface

The screenshot shows the interface with three main panels highlighted by red boxes and numbers:

- 1. The Sidebar Panel (Control Panel):** Contains the 'Upload' section with a 'Data table type' dropdown (set to 'mzMine'), an 'Upload Peak Table (*.csv)' section with a 'Browse...' button, and 'Global parsing settings' including 'Feature ID (auto-generated)', 'm/z column', 'rt column', and 'How to define sample columns?' options.
- 2. The Main Panel (Display Area):** Displays summary metrics: 'TOTAL SAMPLES DETECTED' (9) and 'TOTAL FEATURES (PEAKS)' (6,167). Below this is a 'Dataset: Example' section with a 'Show 6 entries' dropdown and a search bar. A table of data is shown with columns: row ID, row m/z, row retention time, row ion mobility, row ion mobility unit, row CCS, correlation group ID, annotation network number, best ion, auto MS2 verify, identified by n, partners, and neutral M mass. The table shows 6 rows of data. Below the table is a pagination bar showing 'Showing 1 to 6 of 6,167 entries' and a 'Previous' button.
- 3. The Navigation Bar (Tab Menu):** Located at the top, it contains tabs: 'Upload Data', 'Blank Filters', 'MS Filters', 'QC Filters', 'Peak Filters', 'Final Summary', and 'About'.

GUI is organized into 3 main panels:

1. The Sidebar Panel (Control Panel)

- This is the primary area for user interaction. It houses the input controls where users upload data, configure global settings, define parameters, and execute filtering steps.

2. The Main Panel (Display Area)

- This is where the application renders the output and displays interactive data tables, statistical plots, and provides summary metrics.

3. The Navigation Bar (Tab Menu)

- This allows users to navigate linearly through the sequential steps of the data processing pipeline, from data upload to final summary and export.

Contextual assistance is embedded throughout the interface. Users can access detailed parameter descriptions and recommendations by interacting with the tooltip icons.

On the following pages, yellow numbered circles  indicate the recommended operational order within each interface tab in *Control Panel*.

Red numbered squares  describe widgets in *Display Area*.



Implementation Details

All filter groups are applied sequentially in the order shown in the tab panel: Blank Filters → MS Filters → QC Filters → Peak Filters. Changes made on later tabs do not affect the results of previous tabs. To change the order of applied filters (e.g., 'Peak Filters' before 'Blank Filters'), run the application sequentially several times, using the output peak table from one run as the input for the next.

Quick Start

Press button  in “Upload data” tab to upload attached example dataset and get access to the MetaboCensoR functionality. If you use this example, run the Blank filter first to save memory for further steps. You can keep all parameters by default.

Briefly, it is the LC-MS plant profiling dataset, which is described in [Houriet et al.](#) Briefly, methanol extracts from the plant ashwagandha [*Withania somnifera* (L.) Dunal] together with blanks were analyzed on RP column in DDA positive mode on Thermo Q-Exactive Plus Orbitrap. Raw mzML files were then processed in mzMine in default pre-settings for UPLC-DDA to obtain resulted peak table.



Tab “Upload Data”: Control Panel



Upload Data

Blank Filters

MS Filters

QC Filters

Peak Filters

Final Summary

About

Upload

Data table type: ?
mzMine

Upload Peak Table (*.csv)
Browse... No file selected

Global parsing settings

Feature ID (auto-generated):
feature_id

m/z column:
row m/z

rt column:
row retention time

How to define sample columns?
☐ Auto-detect numeric sample columns
☒ By keyword match
☐ Pick columns manually

Sample column keywords (pick/add multiple)
.mzML .mzXML

Clear dataset

Example dataset

TOTAL SAMPLES DETECTED
9

TOTAL FEATURES (PEAKS)
6,167

Dataset: Example

Show 6 entries

Search:

row ID	row m/z	row retention time	row ion mobility	row ion mobility	row CCS	correlation group ID	annotation network number
1							
2							
3							
4							
5							
6							

Showing

next

1

2

3

4

5

1. Pick the right table type and upload
MetaboCensoR supports mzMine, xcms, MS-DIAL, and Default input tables. If the type is wrong, the app throws a parsing error.

2. Confirm the three core columns
Use the generated selectors to check Feature ID (by default auto-generates numeric vector), m/z, and rt. We recommend to start with default column names.

3. Specify sample intensity columns
Use auto-detect, keyword matching, or manual column selection. Keyword matching is the most transparent because you can inspect the exact suffixes/prefixes used in your export.

4. Clear dataset
Reset all calculations.

5. Example dataset shortcut
The Example dataset button loads the built-in example.

7



Tab “Upload Data”: Display Area



[Upload Data](#)[Blank Filters](#)[MS Filters](#)[QC Filters](#)[Peak Filters](#)[Final Summary](#)[About](#)

Upload

Data table type: ?
mzMine

Upload Peak Table (*.csv)
Browse... No file selected

Global parsing settings

Feature ID (auto-generated):
feature_id

m/z column:
row m/z

rt column:
row retention time

How to define sample columns?
☐ Auto-detect numeric sample columns
☒ By keyword match
☐ Pick columns manually

Sample column keywords (pick/add multiple):
.mzML .mzXML

Clear dataset

Example dataset

TOTAL SAMPLES DETECTED
9

1 Statistics of uploaded dataset

Dataset: Example

Show 6 entries Search:

row ID	row m/z	row retention time	row ion mobility	row ion mobility unit	row CCS	correlation group ID	annotation network number
1	108	148.0277472	0.3691534			31	0
2	181	169.0462455	0.3699988			31	
3	175	128.0197664	0.37076914			31	20
4	231	216.951263					
5	43	234.961351	0.3725722			31	14
6	145	244.9475245	0.37269926			31	

Showing 1 to 6 of 6,167 entries Previous 1 2 3 4 5 ... 1,028 Next



Tab “Blank Filters”: Control Panel

The screenshot shows the 'Blank Filters' tab in the MetaboCensus interface. The top navigation bar includes 'Upload Data', 'Blank Filters' (selected), 'MS Filters', 'QC Filters', 'Peak Filters', 'Final Summary', and 'About'. The left sidebar contains the 'Sample Label' and 'Blank / Media Removal' sections. The main area displays the 'Input table' and 'Output table'.

Sample Label Section:

- Label source: ☒ From sample names (by token), ☐ From metadata CSV, ☐ From custom CSV (one column, no header), ☐ Manual editable table
- Token index: 2
- Token separator: -
- ☒ Show labels table

Blank / Media Removal Section:

- Enable Blank filter: ☒
- Choose blank/media group(s): Blank
- Filtering mode: ☒ Blank < Cutoff * Sample Max, ☐ Drop any feature detected in Blank
- Cutoff (fraction of max experimental mean): 0.1
-
- ☒ Show Blank plot panel
-
-

Input table: Shows 5 entries. The first few rows are:

row	108	148.02774
1	108	148.02774
2	181	169.04624
3	175	128.01976
4	231	216.95126
5	43	234.96133

Output table: Shows 5 entries. The first few rows are:

row ID	108	148.02774
1	108	148.02774
2	181	169.04624
4	231	216.95126
5	43	234.96133
6	145	244.94752

Filtering results — Step 1 (Blank Filters):

Metric	Value
Features before	6167
Removed by Blank Filters	3570
Features after	2597

1. Sample Labels — how to specify groups

Choose a label source:

- From sample name:

Define token and index for extracting labels from sample names.

Example — parse from sample names

Sample name: Sample_WT_1_R-C1_1_5037.mzML Peak area

Separator: _

Token index: 2

Resulting label: WT

- From metadata CSV:

Upload a metadata table with column names, choose the sample-name column and the column to use as Label. Rows are matched by sample name. You can also clean sample name to remove file extension, suffixes, etc.

- From custom CSV:

Upload a one-column CSV without header. Label order must match sample column order.

- Manual editable table:

The editable label table is filled from the current token-based labels. Double-click cells in the Label column to manually correct group names.

2. Show resulted Labels in Display Area

3. Set Up Blank Filters

Enable filter.

Specify blank groups.

Choose between filtering modes:

- Blank < Cutoff * Sample Max

- Drop any feature detected in Blank

4. Plot distribution of filtered peaks

5. Apply or Discard Filtering Step



Tab “Blank Filters”: Display Area

MetaboCampor

Upload Data

Blank Filters

MS Filters

QC Filters

Peak Filters

Final Summary

About

Sample Label

Label source:

From sample names (by token)

From metadata CSV

From custom CSV (one column, no header)

Manual editable table

Token index:

2

Token separator:

-

☒ Show labels table

Blank / Media Removal

Enable Blank filter:

Choose blank/media group(s):

Blank

Filtering mode:

Blank < Cutoff * Sample Max

Drop any feature detected in Blank

Cutoff (fraction of max experimental mean):

0.1

Plot values distribution

☒ Show Blank plot panel

Apply

Clear

Input table

Show

5

 entries

Search:

row ID	mz	rt	row ion mobility	row ion mobility unit	row CCS	correlation group ID	annotat network number
1	108	148.0277472	0.3691534			31	0
2	181	169.0462455	0.3699988			31	
3	175	128.0197664	0.37076914			31	20
4	231	216.9512638				31	
5	43	234.961351	0.3725			31	14

Showing 1 to 5 of 6,167 entries

Previous

1

2

3

4

5

...

1,234

Next

Check matching Sample Names and Group Labels:

Show

6

 entries

Sample	Label
Orbi_Blank_A.mzML Peak area	Blank
Orbi_Blank_B.mzML Peak area	Blank
Orbi_Blank_C.mzML Peak area	Blank
Orbi_Sample_A.mzML Peak area	Sample
Orbi_Sample_B.mzML Peak area	Sample
Orbi_Sample_C.mzML Peak area	Sample

Showing 1 to 6 of 9 entries

Previous

1

2

Next

Check Blank Ratio Distribution:

Blank to Sample Ratio Distribution

Count

Ratio: Mean(Blank) / Max(Mean(Sample)) (log10 scale)

Filtered: 57.9%

Output table — after Blank Filters

Show

5

 entries

Search:

row ID	mz	rt	row ion mobility	row ion mobility unit	row CCS	correlation group ID	annotat network number
1	108	148.0277472	0.3691534			31	0
2	181	169.0462455	0.3699988			31	
4	231	216.9512638	0.371403			31	
5	43	234.961351	0.3725			31	14
6	145	244.9475245	0.37269			31	

Showing 1 to 5 of 2,597 entries

Previous

1

2

3

4

5

...

520

Next

Filtering results — Step 1 (Blank Filters)

Metric	Value
Features before	6167
Removed by Blank Filters	3570
Features after	2597



Tab “MS Filters”: Control Panel

The screenshot shows the 'MS Filters' tab in a software interface. The top navigation bar includes 'Upload Data', 'Blank Filters', 'MS Filters' (selected), 'QC Filters', 'Peak Filters', 'Final Summary', and 'About'. The left sidebar contains six filter sections: A) Isotopes & Dimers, B) Adduct Families, C) Neutral Losses, D) In-Source Fragments, E) Mispicked Ions, and F) Saturated Ions. Each section has a title, a toggle for enabling the filter, and various input fields for parameters like correlation, tolerance, and mass. A 'Show' button is present for each section. The main area displays two tables: 'Input table' and 'Output table - after MS'. The 'Input table' shows columns for row ID, mz, rt, row ion mobility, row ion mobility unit, row CCS, correlation group ID, and annotation number. The 'Output table' shows columns for row ID and mz. A 'Download CSV' button is located below the 'Output table'. A yellow box on the right side of the image contains eight numbered steps for using the filters.

A) Isotopes & Dimers

Enable isotope/dimer collapse ☒

Correlation $\geq$ 0.8

rt tolerance (min) 0.01

m/z tolerance type: ☒ Da ☐ ppm

m/z tolerance (Da) 0.005

Max charge (z_max) 3

Max isotope order (n) 1

Dimer-series max n_d 3

B) Adduct Families

Enable adduct-family collapse ☒

Polarity ☒ positive ☐ negative

Adducts CSV (optional)

Browse... Built-in list

Correlation $\geq$ 0.8

rt tolerance (min) 0.005

m/z tolerance type: ☒ Da ☐ ppm

m/z tolerance (Da) 0.005

Minimum neutral mass 50

☒ Strict RT split inside clusters

C) Neutral Losses

Enable neutral loss collapse ☒

D) In-Source Fragments

Enable in-source fragment collapse ☒

E) Mispicked Ions

Enable mispicked ions collapse ☒

F) Saturated Ions

Enable ringing artifact collapse ☒

☒ Show isotope/dimer groups

☐ Show adduct families

☐ Show adduct statistics

☐ Show neutral loss hits

☐ Show fragment clusters

☐ Show mispicked ions groups

☐ Show ringing artifact families

Apply

Clear

Input table

Show 5 entries

Search:

row ID	mz	rt	row ion mobility	row ion mobility unit	row CCS	correlation group ID	annotation number
1	108	148.0277472	0.3691534			31	0
2	181	169.0462455	0.3699988			31	
3	175	128.0197664	0.37076914			31	20
4	151	125.12638	0.371403			31	
5	1351	0.37257					

Output table - after MS

Show 5 entries

row ID	mz
1	108
3	175
7	177
10	171
11	141

Showing 1 to 5 of 4,258 entries

Isotopes table:

Download CSV

Feature	mz
981.5019	5.1
982.5019	5.1
977.9383	9.2
977.4369	9.2
972.4807	9.2
972.9820	9.2
895.9259	4.4
896.4304	4.4
883.8535	4.4
884.3593	4.4

Showing 1 to 5 of 6,167 entries

1. Enable & Set Up Isotopes Filtering

Need parameter help? Click the tooltips or read the [About](#) tab for details. **Pro-Tip:** Start with the default settings.

2. Enable & Set Up Adducts Filtering

Need parameter help? Click the tooltips or read the [About](#) tab for details. **Pro-Tip:** Start with the default settings and check your polarity.

3. Enable & Set Up Neutral Losses Filtering

Need parameter help? Click the tooltips or read the [About](#) tab for details. **Pro-Tip:** Start with the default settings and check your polarity.

4. Enable & Set Up In-Source Fragments Filtering

Need parameter help? Click the tooltips or read the [About](#) tab for details. **Pro-Tip:** Start with the default settings.

4. Enable & Set Up Mispicked Ions Filtering

Need parameter help? Click the tooltips or read the [About](#) tab for details.

6. Enable & Set Up Saturated Ions Filtering

Need parameter help? Click the tooltips or read the [About](#) tab for details.

7. Print & Download annotation tables

8. Apply or Discard Filtering Step



12



Tab “QC Filters”: Control Panel

The screenshot shows the QC Filters control panel with the following components:

- 1. Sample Labels — how to specify groups**: A section for selecting the label source. Options include: Use labels from Blank Tab (selected), From sample names (by token), From metadata CSV, From custom CSV (one column, no header), and Manual editable table. A checkbox for "Show labels table" is checked.
- 2. Show resulted Labels in Display Area**: A section for selecting the label source. Options include: Use labels from Blank Tab (selected), From sample names (by token), From metadata CSV, From custom CSV (one column, no header), and Manual editable table. A checkbox for "Show labels table" is checked.
- 3. Set Up QC Values Filters**: A section for selecting the filter type. Options include: Zeros, Mean, RSD, and Min. A checkbox for "Show QC plot panel" is checked.
- 4. Set Up Zero Values Filters**: A section for selecting the filter type. Options include: Zeros, Mean, RSD, and Min. A checkbox for "Show QC plot panel" is checked.
- 5. Set Up Mean/RSD/Min Values Filters**: A section for selecting the filter type. Options include: Zeros, Mean, RSD, and Min. A checkbox for "Show QC plot panel" is checked.
- 6. Plot distribution of selected values**: A section for selecting the filter type. Options include: Zeros, Mean, RSD, and Min. A checkbox for "Show QC plot panel" is checked.
- 7. Apply or Discard Filtering Step**: A section for selecting the filter type. Options include: Zeros, Mean, RSD, and Min. A checkbox for "Show QC plot panel" is checked.

The panel also includes a "Sample Label" input field, a "QC value filters" section with "Enable filters" and "Evaluation mode" (Groupwise, EVERY, POOLED), a "Group(s) to evaluate" field, a "Zeros metric" section (Count, Percent), a "Zeros ≤ threshold (count or %)" field, a "Mean ≥ threshold" field, an "RSD ≤ threshold" field, a "Min ≥ threshold" field, a "Plot values distribution" checkbox, and "Apply" and "Clear" buttons.

The "Output table" shows the results of the filtering process. The table has columns for row ID, mz, and various metrics. The "Output table — after QC" shows the results after applying the filters. The table has columns for row ID, mz, and various metrics. The "Filtering results — Step 3 (QC Filters)" table shows the metrics and values for each step of the filtering process.

Metric	Value
Features before	2597
Removed by Zeros	5
Removed by Mean	1453
Removed by RSD	110
Removed by Min	595
Total removed	1485
Features after	1112

1. Sample Labels — how to specify groups

Choose a label source: from *Blank Filters* tab, sample names, metadata CSV, custom CSV (one column, no header), or by manual editing.

See details for extracting Labels on the Page **Tab “Blank Filters”: Control Panel**.

2. Show resulted Labels in Display Area

3. Set Up QC Values Filters

Select desired values for filtering:
Zeros, Mean, RSD, or Min.

Choose between filtering modes:

- Groupwise (ANY selected group passes)
- Groupwise (EVERY selected group passes)
- POOLED across selected groups

Specify groups used for evaluation.

4. Set Up Zero Values Filters

Specify Zero metric (counts or %) and threshold.

5. Set Up Mean/RSD/Min Values Filters

Specify thresholds for corresponding value.

6. Plot distribution of selected values

7. Apply or Discard Filtering Step



14



Tab “Peak Filters”: Control Panel

The screenshot shows the 'Peak Filters' control panel with the following components:

- Navigation Bar:** Upload Data, Blank Filters, MS Filters, QC Filters, **Peak Filters**, Final Summary, About.
- Peak Filters Panel:**
 - Enable filters:** m/z, rt, RMD, AMD (all checked).
 - Target peak list:** A text input field.
 - m/z filters:** m/z ≥ 200, m/z ≤ (empty).
 - rt filters:** rt ≥ (min) 1, rt ≤ (min) (empty).
 - Action for matched target peaks:** Remove matched peaks (selected), Keep only matched peaks.
 - m/z tolerance type:** Da (selected), ppm.
 - m/z tolerance (Da):** 0.005.
 - rt tolerance (min):** 0.05.
 - Type target peaks:** 506.528, 11.6, 970.611, 9.7.
 - Or upload target list (.csv):** Browse... No file selected.
 - CSV note:** CSV should contain mz and rt columns. If column names are not detected, the first two columns are treated as mz and rt.
 - Plot values distribution:** Show Peak plot panel (checked).
 - Buttons:** Apply, Clear.
- Input table:** A table with columns: row ID, mz, rt, row ion mobility, row ion mobility unit, correlation group ID, annotation number. It shows 5 entries.
- Matched Peaks:** A table with columns: target_id, target_mz, target_rt, Feature, mz, rt, delta_mz, delta_mz_ppm. It shows 2 entries.
- Scatter Plot:** A plot of m/z vs. RT (min) with a red dashed line at m/z 252. A label 'Filtered m/z: 252' is present.
- Output table – after Peak Filters:** A table with columns: row ID, mz, rt, row ion mobility. It shows 5 entries.
- Filtering results – Step 4 (Peak Filters):** A table with columns: Metric, Value.

Metric	Value
Features before	2597
Removed by m/z	252
Removed by rt	555
Removed by RMD	0
Removed by AMD	0
Matched target peaks	2
Removed by target m/z + rt	2
Total removed	656
Features after	1941

Note: a feature may be counted in multiple filter categories.

1. Set Up Peak Values Filters

Select desired values for filtering:

m/z, rt, mass defect (RMD, or AMD), or a list of target peaks to remove.

2. Set Up m/z/rt/AMD/RMD Values Filters

Specify threshold ranges for corresponding value.

3. Set Up Peaks to Remove/Keep

Specify action as keep or remove; provide target m/z and RT pairs. Enter comma-separated values manually, one pair per line, or upload a CSV file. Specify m/z and RT tolerances for matching.

4. Plot distribution of selected values

5. Apply or Discard Filtering Step



Tab “Peak Filters”: Display Area

MC

ADAMC

Upload DataBlank FiltersMS FiltersQC FiltersPeak FiltersFinal SummaryAbout

Peak Filters

Enable filters: 1

m/zrtRMDAMD

Target peak list

☒ m/z ≥

☐ m/z ≤

200

☒ rt ≥ (min)

☐ rt ≤ (min)

1

Action for matched target peaks:

☒ Remove matched peaks

☐ Keep only matched peaks

m/z tolerance type:

☒ Da

☐ ppm

m/z tolerance (Da)

0.005

rt tolerance (min)

0.05

Type target peaks:

506.528, 11.6

970.611, 9.7

Or upload target list (.csv)

Browse...

No file selected

CSV should contain mz and rt columns. If column names are not detected, the first two columns are treated as mz and rt.

Plot values distribution

☒ Show Peak plot panel

Apply

Clear

Input table

Show 5 entries

Search:

row ID	mz	rt	row ion mobility	row CCS	correlation group ID	annotation number
1	108	148.0277472	0.3691534		31	0
2	181	169.0462455	0.3699988		31	
4	231	216.951263			31	
5	43	234.96135			31	14
6	145	244.9475245	0.37269926		31	

Showing 1 to 5 of 2,597 entries

Previous12345...520Next

Matched Peaks

Show 5 entries

Search:

target_id	target_mz	target_rt	mz	rt	delta_mz	delta_mz_ppm	
typed_1	506.5280	11.600	24	506.5288	11.598	0.0008	1.5
ty							0.4
1.4							1.2

Showing 1 to 5 of 2,597 entries

Previous12345...520Next

Check values map:

m/z vs RT

1500

1000

500

0

3

6

9

12

RT (min)

Filtered m/z: 252

Filtered rt: 555

Total Filtered: 654

Output table – after Peak Filters

Show 5 entries

Search:

row ID	mz	rt	row ion mobility	row CCS	correlation group ID	annotation number
818	2815	336.074296	1.0103563		58	489
819	2820	314.0921816	1.0123757		58	489
820	2814	284.134450				490
821	2819	289.09139				490
822	2824	352.0481904				489

Showing 1 to 5 of 1,941 entries

Previous12345...389Next

Filtering results – Step 4 (Peak Filters)

Metric	Value
Features before	2597
Removed by m/z	252
Removed by rt	555
Removed by RMD	0
Removed by AMD	0
Matched target peaks	2
Removed by target m/z + rt	2
Total removed	656
Features after	1941

16



Tab “Final Summary”: Control Panel



Upload Data

Blank Filters

MS Filters

QC Filters

Peak Filters

Final Summary

About

Compile final output

Compile output ?

Final peak table (.csv)

MGF Filtering

Filter MGF ☒ ?

Upload MGF file

Browse...

orbi_jimn_gms.mgf

Upload complete

Select ID column in Table:

row ID

MGF Feature ID field:

acquisitionNum

Filter MGF

Clear MGF

Success! Matched 321 spectra from your filtered table.

Download Filtered MGF

Filtering summary

- Dataset: Example
- Step 1 Blank Filters: APPLIED
- Step 2 MS Filters: APPLIED

1. Compile final filtered table

Clicking this button enables the table download and MGF filtering.

2. Download final table

3. Enable MGF filtering & upload .mgf file

4. Match peak ID between table and .mgf file

The feature identifier in the final peak table must correspond to the chosen spectrum identifier field inside the uploaded MGF. We recommend using the default columns unless you have specific requirements.

5. Filter or clear uploaded .mgf file

6. Download filtered .mgf file

Showing 1 to 5 of 1,083 entries

Previous

1

2

3

4

5

...

217

Next

	annotat	
	network	
	number	
819		489
821		490
822		489
829		495
830		495



Tab “Final Summary”: Display Area

[Upload Data](#)[Blank Filters](#)[MS Filters](#)[QC Filters](#)[Peak Filters](#)[Final Summary](#)[About](#)

Compile final output

[Compile output](#)[Final peak table \(.csv\)](#)

MGF Filtering

Filter MGF ☒

Upload MGF file

[Browse...](#)

orbi_iimn_gnps.mgf

Upload complete

Select ID column in Table:

row ID

MGF Feature ID field:

acquisitionNum

[Filter MGF](#)[Clear MGF](#)

Success! Matched 321 spectra from your filtered table.

[Download Filtered MGF](#)

Filtering summary

- Dataset: Example
- Step 1 Blank Filters: APPLIED
- Step 2 MS Filters: APPLIED
- Step 3 QC Filters: not applied
- Step 4 Peak Filters: APPLIED
- Raw features: 6167
- Filtered features: 5084 (82.44%)
- Final features: 1083

Final table

Show 5 entries

Search:

row ID	row m/z	row retention time	row ion mobility	row ion mobility unit	row CCS	correlation group ID	annotation network number
819	2820	314.0921816	1.0127777			58	489
821	2819	289.091359	1.016			58	490
822	2824	352.0481904	1.0176588			58	489
829	2893	275.1104506	1.0585195			60	495
830	2862	253.1291676	1.0596385			60	495

Showing 1 to 5 of 1,083 entries

[Previous](#)[1](#)[2](#)[3](#)[4](#)[5](#)[...](#)[217](#)[Next](#)

3 Matches in filtered .mgf



Tab “About” Overview

How to:

Select a section:

- Quick Start
- Upload Data
- Blank Filters
- MS Filters
- QC Filters
- Peak Filters
- Final Summary
- Troubleshooting
- Project Details

Troubleshooting

- **No sample columns:** your sample keywords don't match sample column names.
- **Parsing error:** check selected Data table type and compatibility with examples of accessible Peak Table data formats in the [GitHub](#) repository.
- **Disconnected from the server:** could be caused by memory limit, try to compile and run it locally. Instructions in the [GitHub](#) repository.

Navigate through the sections to get more information

by Plyushchenko I.V. | 2026 | [GitHub](#) | [Cite](#)