



MetaboCensoR Tutorial

2026-06-03

Tutorial Edition 1.1

Software Release 1.1.x



Contents

General Information.....	3
Running & Installation.....	4
Input Requirements.....	4
Overview of Interface	5
Implementation Details.....	6
Quick Start	6
Tab “Upload Data”: Control Panel.....	7
Tab “Upload Data”: Display Area.....	8
Tab “Blank Filters”: Control Panel	9
Tab “Blank Filters”: Display Area	10
Tab “MS Filters”: Control Panel.....	11
Tab “MS Filters”: Display Area	12
Tab “QC Filters”: Control Panel	13
Tab “QC Filters”: Display Area	14
Tab “Peak Filters”: Control Panel.....	15
Tab “Peak Filters”: Display Area	16
Tab “Final Summary”: Control Panel	17
Tab “Final Summary”: Display Area	18
Tab “About” Overview	19



General Information

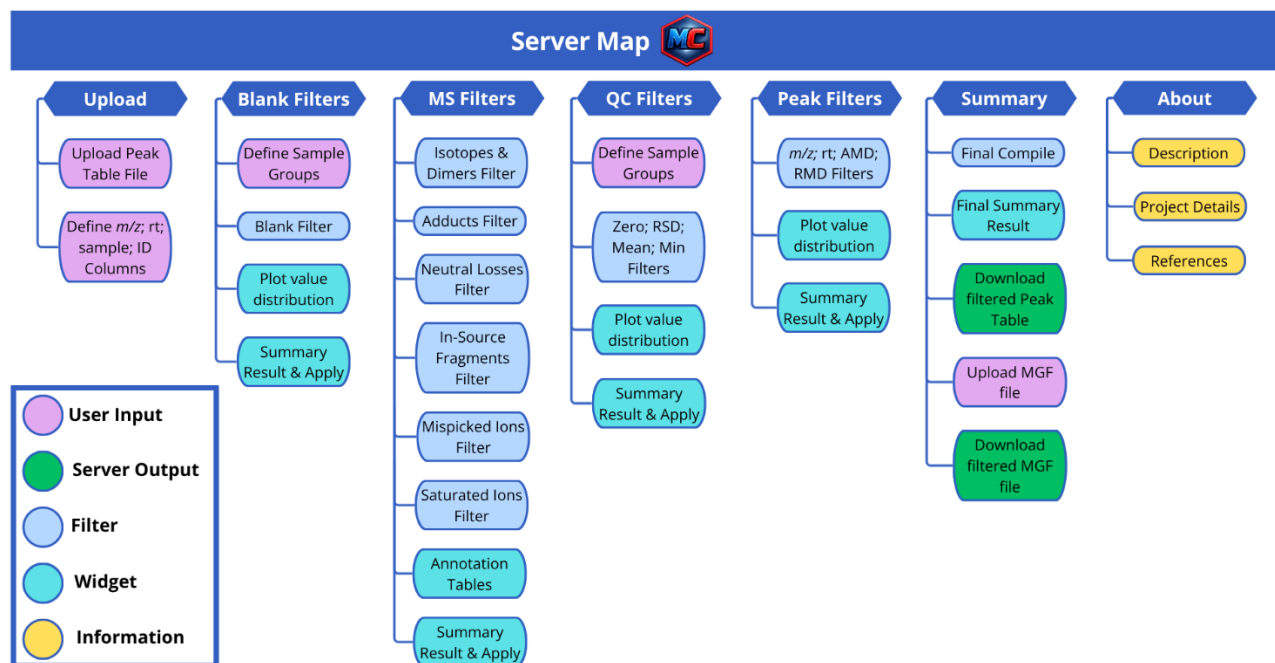
This tutorial aims at guiding the user to process data LC-MS data through Shiny GUI interface of MetaboCensoR. MetaboCensoR is an input-versatile Shiny application and local R package for analyte-centric peak table filtering. The workflow integrates four complementary modules for blank filtering, redundant ion-species filtering, quality-control filtering, and peak-based filtering. MetaboCensoR also provides interactive threshold optimization, exportable annotation tables, and synchronized filtering of associated .mgf files.

An overview of MetaboCensoR is presented below. The workflow is organized into six functional tabs that collectively support peak table curation, parameter optimization, and export of filtered results. The output of the app is a filtered peak table and, optionally, an MGF spectra file. Only the peak table upload is required for full workflow execution and final output generation. All other tabs are optional.

License: GPL-3.0

Contact: plyushchenko.ivan@gmail.com

Project Repository: github.com/plyush1993/MetaboCensoR





Running & Installation

Shiny Web Deployment is available by the link:

plyush1993.shinyapps.io/metabocensor/

Local installation of R package:

```
if (!require("BiocManager", quietly = TRUE)) {  
  install.packages("BiocManager")  
}  
if (!require("remotes", quietly = TRUE)) {  
  install.packages("remotes")  
}  
remotes::install_github("plyush1993/MetaboCensor", INSTALL_opts = "--no-multiarch")
```

or

```
if (!requireNamespace("pak", quietly = TRUE)) install.packages("pak")  
pak::pak("plyush1993/MetaboCensor")
```

Run Locally:

```
MetaboCensor::run_metabocensor()
```

Input Requirements

The application accepts peak tables in **.csv** format and MS1/MS2 mass spectral data in **.mgf** format, both of which can be directly exported from dedicated metabolomics processing software (xcms, MZmine, MS-DIAL). Make sure that your input peak table is aligned with data formats. Examples of acceptable formats are available here: [Input Examples](#).

Important note

Although the web version has been successfully tested on datasets with over 10,000 peaks, the Shiny server has limited memory capacity. Computationally intensive steps (like MS filters) or large MGF files (> 3 MB) may cause the app to stack or disconnect. For heavy processing, you should run the application locally.



Overview of Interface

The screenshot displays the MC interface with three main panels highlighted by red boxes and numbers:

- 1. Sidebar Panel (Control Panel):** Contains an 'Upload' section with a 'Data table type' dropdown (set to 'mzMine') and an 'Upload Peak Table (*.csv)' section with a 'Browse...' button. Below this is a 'Global parsing settings' section with dropdowns for 'Feature ID (auto-generated)', 'm/z column', and 'rt column'. At the bottom, there are radio buttons for 'How to define sample columns?': 'Auto-detect numeric sample columns', 'By keyword match' (selected), and 'Pick columns manually'.
- 2. Main Panel (Display Area):** Displays a dataset table titled 'Dataset: Example'. The table has columns: row ID, row m/z, row retention time, row ion mobility, row ion mobility unit, row CCS, correlation group ID, annotation network number, best ion, auto MS2 verify, identified by n, partners, and neutral M mass. The table shows 6 entries. Above the table, it says 'TOTAL SAMPLES DETECTED: 9' and 'TOTAL FEATURES (PEAKS): 6,167'. Below the table, it says 'Showing 1 to 6 of 6,167 entries' and has pagination controls.
- 3. Navigation Bar (Tab Menu):** Located at the top, it contains tabs: 'Upload Data', 'Blank Filters', 'MS Filters', 'QC Filters', 'Peak Filters', 'Final Summary', and 'About'.

GUI is organized into 3 main panels:

1. The Sidebar Panel (Control Panel)

- This is the primary area for user interaction. It houses the input controls where users upload data, configure global settings, define parameters, and execute filtering steps.

2. The Main Panel (Display Area)

- This is where the application renders the output and displays interactive data tables, statistical plots, and provides summary metrics.

3. The Navigation Bar (Tab Menu)

- This allows users to navigate linearly through the sequential steps of the data processing pipeline, from data upload to final summary and export.

Contextual assistance is embedded throughout the interface. Users can access detailed parameter descriptions and recommendations by interacting with the tooltip icons.

On the following pages, yellow numbered circles  indicate the recommended operational order within each interface tab in *Control Panel*.

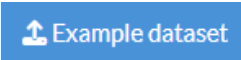
Red numbered squares  describe widgets in *Display Area*.



Implementation Details

All filter groups are applied sequentially in the order shown in the tab panel: Blank Filters → MS Filters → QC Filters → Peak Filters. Changes made on later tabs do not affect the results of previous tabs. To change the order of applied filters (e.g., 'Peak Filters' before 'Blank Filters'), run the application sequentially several times, using the output peak table from one run as the input for the next.

Quick Start

Press button  in “Upload data” tab to upload attached example dataset and get access to the MetaboCensoR functionality. If you use this example, run the Blank filter first to save memory for further steps. You can keep all parameters by default.

Briefly, it is the LC-MS plant profiling dataset, which is described in [Houriet et al.](#) Briefly, methanol extracts from the plant ashwagandha [*Withania somnifera* (L.) Dunal] together with blanks were analyzed on RP column in DDA positive mode on Thermo Q-Exactive Plus Orbitrap. Raw mzML files were then processed in mzMine in default pre-settings for UPLC-DDA to obtain resulted peak table.



Tab “Upload Data”: Control Panel



Upload Data

Blank Filters

MS Filters

QC Filters

Peak Filters

Final Summary

About

Upload

Data table type: ?
mzMine

Upload Peak Table (*.csv)
Browse... No file selected

Global parsing settings

Feature ID (auto-generated):
feature_id

m/z column:
row m/z

rt column:
row retention time

How to define sample columns?
☐ Auto-detect numeric sample columns
☒ By keyword match
☐ Pick columns manually

Sample column keywords (pick/add multiple)
.mzML .mzXML

Clear dataset

Example dataset

?

TOTAL SAMPLES DETECTED
9

TOTAL FEATURES (PEAKS)
6,167

Dataset: Example

Show 6 entries

Search:

row ID	row m/z	row retention time	row ion mobility	row ion mobility	row CCS	correlation group ID	annotation network number
1							
2							
3							
4							
5							
6							

Showing

next

1

Pick the right table type and upload

MetaboCensoR supports mzMine, xcms, MS-DIAL, and Default input tables. If the type is wrong, the app throws a parsing error.

2

Confirm the three core columns

Use the generated selectors to check Feature ID (by default auto-generates numeric vector), m/z, and rt. We recommend to start with default column names.

3

Specify sample intensity columns

Use auto-detect, keyword matching, or manual column selection. Keyword matching is the most transparent because you can inspect the exact suffixes/prefixes used in your export.

4

Example dataset shortcut

The Example dataset button loads the built-in example.

by Plyushchenko I.V. | 2026 | GitHub | Cite



Tab “Upload Data”: Display Area



[Upload Data](#)[Blank Filters](#)[MS Filters](#)[QC Filters](#)[Peak Filters](#)[Final Summary](#)[About](#)

Upload

Data table type: ?
mzMine

Upload Peak Table (*.csv)
Browse... No file selected

Global parsing settings

Feature ID (auto-generated):
feature_id

m/z column:
row m/z

rt column:
row retention time

How to define sample columns?
☐ Auto-detect numeric sample columns
☒ By keyword match
☐ Pick columns manually

Sample column keywords (pick/add multiple):
.mzML .mzXML

Clear dataset

Example dataset

TOTAL SAMPLES DETECTED
9

1 Statistics of uploaded dataset

Dataset: Example

Show 6 entries Search:

row ID	row m/z	row retention time	row ion mobility	row ion mobility unit	row CCS	correlation group ID	annotation network number
1	108	148.0277472	0.3691534			31	0
2	181	169.0462455	0.3699988			31	
3	175	128.0197664	0.37076914			31	20
4	231	216.951263					
5	43	234.961351	0.3725722			31	14
6	145	244.9475245	0.37269926			31	

Showing 1 to 6 of 6,167 entries Previous 1 2 3 4 5 ... 1,028 Next



Tab “Blank Filters”: Control Panel



Sample Label

Label source: ?
☒ From sample names
☐ From custom CSV (one column, no header)

Token index
2

Token separator
-

☒ Show labels table

Blank / Media Removal

Enable Blank filter ☒

Choose blank/media group(s):
Pick one or more (e.g., Media, Blank)...

Filtering mode:
☒ Blank < Cutoff * Sample Max
☐ Drop any feature detected in Blank

Cutoff (fraction of max experimental mean):
0.1

Plot values distribution
☒ Show Blank plot panel

Apply

Clear

Input table

Show 5 entries

Search:

row ID	mz	rt	row ion mobility	row ion mobility unit	row CCS	correlation group ID	annotation network number
1	108	148.0277472	0.3691534			31	0
2	181	169.0462455	0.3699988			31	
3	175	128.0197664	0.37076914			31	20
4	231	216.9512638	0.371403			31	
5	43	234.961351	0.3725722			31	14

1. Sample Labels — how to specify groups

Choose a label source: from sample names or custom CSV (one column, no header).

Define token and index for extracting labels from sample names.

Example — parse from sample names

Sample name: Sample_WT_1_R-C1_1_5037.mzML Peak area

Separator: _

Token index: 2

Resulting label: WT

2. Show resulted Labels in Display Area

3. Set Up Blank Filters

Enable filter.

Specify blank groups.

Choose between filtering modes:

- Blank < Cutoff * Sample Max
- Drop any feature detected in Blank

4. Plot distribution of filtered peaks

5. Apply or Discard Filtering Step

1							
2							
3	175	128.0197664	0.37076914			31	20
4	231	216.9512638	0.371403			31	
5	43	234.961351	0.3725722			31	14



Tab “Blank Filters”: Display Area

Upload Data

Blank Filters

MS Filters

QC Filters

Peak Filters

Final Summary

About

Sample Label

Label source:

☒ From sample names

☐ From custom CSV (one column, no header)

Token index

2

Token separator

-

☒ Show labels table

Blank / Media Removal

Enable Blank filter

Choose blank/media group(s):

Blank

Filtering mode:

☒ Blank < Cutoff * Sample Max

☐ Drop any feature detected in Blank

Cutoff (fraction of max experimental mean):

0.1

Plot values distribution

☒ Show Blank plot panel

Apply

Clear

Input table

Show 5 entries

Search:

row ID	mz	rt	row ion mobility	row ion mobility unit	row CCS	correlation group ID	annotatic network number
1	108	148.0277472	0.3691534			31	0
2	181	169.0462455	0.3699988			31	
3	175	128.0197664	0.37076914			31	20
4	231	216.9512638				31	
5	43	234.961351				31	14

Showing 1 to 5 of 6,167 entries

Previous12345...1,234Next

Check matching Sample Names and Group Labels:

Show 6 entries

Search:

Sample	Label
Orbi_Blank_A.mzML Peak area	Blank
Orbi_Blank_B.mzML Peak area	Blank
Orbi_Blank_C.mzML Peak area	Blank
Orbi_Sample_A.mzML Peak area	Sample
Orbi_Sample_B.mzML Peak area	Sample
Orbi_Sample_C.mzML Peak area	Sample

Showing 1 to 6 of 9 entries

Previous12Next

Check Blank Ratio Distribution:

Blank to Sample Ratio Distribution

Count

Ratio: Mean(Blank) / Max(Mean(Sample)), (log10 scale)

Filtered: 57.9%

5e+041.5

Output table — after Blank Filters

Show 5 entries

Search:

row ID	mz	rt	row ion mobility	row ion mobility unit	row CCS	correlation group ID	annotatic network number
1	108	148.0277472	0.3691534			31	0
2	181	169.0462455	0.3699988			31	
4	231	216.9512638	0.37076914			31	
5	43	234.961351	0.37269926			31	14
6	145	244.9475245	0.37269926			31	

Showing 1 to 5 of 2,597 entries

Previous12345...520Next

Filtering results — Step 1 (Blank filters)

Metric

Value

Features before

6167

Removed by Blank Filters

3570

Features after

2597



Tab “MS Filters”: Control Panel

Upload Data | Blank Filters | **MS Filters** | QC Filters | Peak Filters | Final Summary | About

A) Isotopes & Dimers

Enable isotope/dimer collapse ☒

Correlation $\geq$ 0.8

rt tolerance (min) 0.01

m/z tolerance type: ☒ Da ☐ ppm

m/z tolerance (Da) 0.005

Max charge (z_max) 3

Max isotope order (n) 1

Dimer-series max n_d 3

B) Adduct Families

Enable adduct-family collapse ☒

Polarity ☒ positive ☐ negative

Adducts CSV (optional)

Browse... Built-in list

Correlation $\geq$ 0.8

rt tolerance (min) 0.005

m/z tolerance type: ☒ Da ☐ ppm

m/z tolerance (Da) 0.005

Minimum neutral mass 50

☒ Strict RT split inside clusters

C) Neutral Losses

Enable neutral loss collapse ☐

D) In-Source Fragments

Enable in-source fragment collapse ☐

E) Mispicked Ions

Enable mispicked ions collapse ☐

F) Saturated Ions

Enable ringing artifact collapse ☐

☒ Show isotope/dimer groups

☐ Show adduct families

☐ Show adduct statistics

☐ Show neutral loss hits

☐ Show fragment clusters

☐ Show mispicked ions groups

☐ Show ringing artifact families

Apply

Clear

Input table

Show 5 entries

Search:

row ID	mz	rt	row ion mobility	row ion unit	row CCS	correlation group ID	annotation number
1	108	148.0277472	0.3691534			31	0
2	181	169.0462455	0.3699988			31	
3	175	128.0197664	0.37076914			31	20
4	151	128.0197664	0.371403			31	
5	1351	0.37257					

Output table — after MS

Show 5 entries

row ID	mz
1	108
3	175
7	177
10	171
11	141

Showing 1 to 5 of 4,258 entries

Isotopes table:

Download CSV

Feature	mz
981.5019	5.5
982.5019	5.5
977.9383	9.2
977.4369	9.2
972.4807	9.2
972.9820	9.2
895.9259	4.4
896.4304	4.4
883.8535	4.4
884.3593	4.4

Showing 1 to 5 of 6,167 entries

1. Enable & Set Up Isotopes Filtering

Need parameter help? Click the tooltips or read the [About](#) tab for details. **Pro-Tip:** Start with the default settings.

2. Enable & Set Up Adducts Filtering

Need parameter help? Click the tooltips or read the [About](#) tab for details. **Pro-Tip:** Start with the default settings and check your polarity.

3. Enable & Set Up Neutral Losses Filtering

Need parameter help? Click the tooltips or read the [About](#) tab for details. **Pro-Tip:** Start with the default settings and check your polarity.

4. Enable & Set Up In-Source Fragments Filtering

Need parameter help? Click the tooltips or read the [About](#) tab for details. **Pro-Tip:** Start with the default settings.

4. Enable & Set Up Mispicked Ions Filtering

Need parameter help? Click the tooltips or read the [About](#) tab for details.

6. Enable & Set Up Saturated Ions Filtering

Need parameter help? Click the tooltips or read the [About](#) tab for details.

7. Print & Download annotation tables

8. Apply or Discard Filtering Step



12



Tab “QC Filters”: Control Panel

The screenshot shows the QC Filters Control Panel with the following components and numbered callouts:

- 1. Sample Labels — how to specify groups**: A section for specifying sample labels, including a search bar and a table of sample names.
- 2. Show resulted Labels in Display Area**: A checkbox to show the resulting labels in the display area.
- 3. Set Up QC Values Filters**: A section for setting up QC value filters, including a table of filter values.
- 4. Set Up Zero Values Filters**: A section for setting up zero value filters, including a table of filter values.
- 5. Set Up Mean/RSD/Min Values Filters**: A section for setting up mean, RSD, and min value filters, including a table of filter values.
- 6. Plot distribution of selected values**: A section for plotting the distribution of selected values, including a histogram.
- 7. Apply or Discard Filtering Step**: A section for applying or discarding the filtering step, including an 'Apply' button.

1. Sample Labels — how to specify groups

Choose a label source: from *Blank Filters* tab, sample names or custom CSV (one column, no header).

See details for extracting Labels on the Page **Tab “Blank Filters”: Control Panel**.

2. Show resulted Labels in Display Area

Select desired values for filtering:
Zeros, Mean, RSD, or Min.

Choose between filtering modes:

- Groupwise (ANY selected group passes)
- Groupwise (EVERY selected group passes)
- POOLED across selected groups

Specify groups used for evaluation.

4. Set Up Zero Values Filters

Specify Zero metric (counts or %) and threshold.

5. Set Up Mean/RSD/Min Values Filters

Specify thresholds for corresponding value.

6. Plot distribution of selected values

7. Apply or Discard Filtering Step

Output table — after QC Filters

Show 5 entries

row ID	mz	rt	row ion mobility	row ion mobility unit	row CCS	correlation group ID	annotatic network number
1	108	148.0277472	0.3691534			31	0
2	181	169.0462455	0.3699988			31	
3	175	128.0197664	0.37076914			31	20
4	231	216.9512638	0.371403			31	
5	43	234.961351	0.3725722			31	14

Showing 1 to 5 of 1,881 entries

Filtering results — Step 3 (QC Filters)

Metric	Value
Features before	6129
Removed by Zeros	980
Removed by Mean	4083
Removed by RSD	2254



Tab “QC Filters”: Display Area

Sample Label

Label source:

Use labels from Blank step

From sample names

From custom CSV (one column, no header)

Show labels table

QC value filters

Enable filters:

Zeros

Mean

RSD

Min

Evaluation mode:

Groupwise (ANY selected group passes)

Groupwise (EVERY selected group passes)

POOLED across selected groups

Group(s) to evaluate:

Sample

Zeros metric

Count

Percent

Zeros $\leq$ threshold (count or %)

0

Mean $\pm$ threshold

5000

RSD $\leq$ threshold

30

Min $\pm$ threshold

1000

Plot values distribution

Show QC plot panel

Apply

Clear

Input table

Show 5 entries

Search:

row ID	mz	rt	row ion mobility unit	row CCS	correlation group ID	annotatic network number
1	108	148.0277472	0.3691534		31	0
2	181	169.0462			31	
3	175	128.0197			31	20
4	231	216.9512638	0.371403		31	
5	43	234.961351	0.3725722		31	14

Showing 1 to 5 of 6,129 entries

Previous 1 2 3 4 5 ... 1,226 Next

Sample-Label table

Show 6 entries

Search:

Sample	Label
Orbi_Blank_A,mzML Peak area	Blank
Orbi_Blank_B,mzML Peak area	Blank
Orbi_Blank_C,mzML Peak area	Blank
Orbi_S	
Orbi_S	
Orbi_Sample_C,mzML Peak area	Sample

Showing 1 to 6 of 9 entries

Previous 1 2 Next

Distribution plots

Check Values Distribution:

ZEROS [Zero values (count)] - Groupwise (ANY) | Filtered: 980 of 6129 (16%)

MEAN [Mean intensity (log10)] - Groupwise (ANY) | Filtered: 4083 of 6129 (66.6%)

MIN [Minimum intensity (log10)] - Groupwise (ANY) | Filtered: 2742 of 6129 (44.7%)

RSD [RSD (%)] - Groupwise (ANY) | Filtered: 1618 of 5493 (29.5%)

Output table — after QC Filters

Show 5 entries

Search:

row ID	mz	rt	row ion mobility unit	row CCS	correlation group ID	annotatic network number
1	108	148.0277472	0.3691534		31	0
2	181	169.0462	0.2698889		31	
3	175	128.0197			31	20
4	231	216.95			31	
5	43	234.961351	0.3725722		31	14

Showing 1 to 5 of 1,881 entries

Previous 1 2 3 4 5 ... 377 Next

Summary

Filtering results — Step 3 (QC Filters)

Metric	Value
Features before	6129
Removed by Zeros	980
Removed by Mean	4083
Removed by RSD	2254
Removed by Min	2742
Total removed	4248
Features after	1881

Note: a feature may be counted in multiple filter categories.

14



Tab “Peak Filters”: Control Panel

Peak Filters

Enable filters: **m/z** **rt** RMD AMD

☒ m/z ≥ 200 ☐ m/z ≤

☒ rt ≥ 1 ☐ rt ≤

Plot values distribution
☒ Show Peak plot panel

Apply

Clear

Input table

Show 5 entries Search:

row ID	m/z	rt	row ion mobility	row ion mobility unit	row CCS	correlation group ID	annotation network number
1	108	148.0277472	0.3691534			31	0
2	169	0.462455	0.3699988			31	
3	216	9512638	0.371403			31	
4	234	961351	0.3725722			31	14
5	6	145	244.9475245	0.37269926		31	

Showing 1 to 5 of 1,068 entries Previous 1 2 3 4 5 ... 214 Next

m/z vs RT

Filtered m/z: 157
Filtered rt: 316
Total Filtered: 369

1. Set Up Peak Values Filters

Select desired values for filtering:
m/z, *rt*, or mass defect (RMD, or AMD).

2. Set Up *m/z*/*rt*/AMD/RMD Values Filters

Specify threshold ranges for corresponding value.

3. Plot distribution of selected values

4. Apply or Discard Filtering Step

Output table — after Peak Filters

Show 5 entries Search:

row ID	m/z	rt	row ion mobility	row ion mobility unit	row CCS	correlation group ID	annotation network number
818	2815	336.074296	1.0103563			58	489
819	2820	314.0921816	1.0123757			58	489
825	2879	291.0843679	1.0544833			60	495
829	2893	275.1104506	1.0585195			60	495
830	2862	253.1291676	1.0596385			60	495

Showing 1 to 5 of 699 entries Previous 1 2 3 4 5 ... 140 Next

Filtering results — Step 4 (Peak Filters)

Metric	Value
Features before	1068
Removed by m/z	157
Removed by rt	316
Removed by RMD	0
Removed by AMD	0
Total removed	369



Tab “Peak Filters”: Display Area

MC

MassCloudR

Upload Data

Blank Filters

MS Filters

QC Filters

Peak Filters

Final Summary

About

Peak Filters

Enable filters: 1

m/z

rt

RMD

AMD

☒ m/z ≥

☐ m/z ≤

200

☒ rt ≥ (min)

☐ rt ≤ (min)

1

Plot values distribution

☒ Show Peak plot panel

Apply

Clear

Input table

Showing 1 to 5 of 1,068 entries

1

2

3

4

5

...

214

Next

Check values map:

m/z vs RT

1500

1000

500

0

0

1000

1500

Filtered m/z: 157

Filtered rt: 316

Total Filtered: 369

Distribution plot

Output table – after Peak Filters

Showing 1 to 5 of 699 entries

1

2

3

4

5

...

140

Next

Filtering results – Step 4 (Peak Filters)

Metric

Value

Features before

1068

Removed by m/z

157

Removed by rt

316

Removed by RMD

0

Removed by AMD

0

Total removed

369

Features after

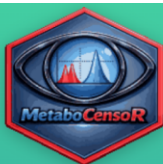
699

Summary

Note: a feature may be counted in multiple filter categories.



Tab “Final Summary”: Control Panel

[Upload Data](#)[Blank Filters](#)[MS Filters](#)[QC Filters](#)[Peak Filters](#)[Final Summary](#)[About](#)

Compile final output

[Compile output](#)[Final peak table \(.csv\)](#)

MGF Filtering

Filter MGF ☒

Upload MGF file

[Browse...](#)

mgf_gnps_test.mgf

Upload complete

Select ID column in Table:

row ID

MGF Feature ID field:

acquisitionNum

[Filter MGF](#)[Clear MGF](#)

Success! Matched 2 spectra from your filtered table.

[Download Filtered MGF](#)

Filtering summary

- Dataset: Example
- Step 1 Blank Filters: APPLIED
- Step 2 MS Filters: APPLIED

1. Compile final filtered table

Clicking this button enables the table download and MGF filtering.

2. Download final table

3. Enable MGF filtering & upload .mgf file

4. Match peak ID between table and .mgf file

The feature identifier in the final peak table must correspond to the chosen spectrum identifier field inside the uploaded MGF. If no IDs match, the app reports: “No matching spectra were found in the MGF file”.

We recommend using the default columns unless you have specific requirements.

5. Filter or clear uploaded .mgf file

6. Download filtered .mgf file



Tab “Final Summary”: Display Area



[Upload Data](#)[Blank Filters](#)[MS Filters](#)[QC Filters](#)[Peak Filters](#)[Final Summary](#)[About](#)

Compile final output

Compile output ?

Final peak table (.csv)

MGF Filtering

Filter MGF ☒ ?

Upload MGF file

Browse...

mgf_gnps_toy.mgf

Upload complete

Select ID column in Table:

row ID

MGF Feature ID field:

acquisitionNum

Filter MGF

Clear MGF

Success! Matched 2 spectra from your filtered table.

Download Filtered MGF

Filtering summary

- Dataset: Example
- Step 1 Blank Filters: APPLIED
- Step 2 MS Filters: APPLIED
- Step 3 QC Filters: not applied
- Step 4 Peak Filters: not applied
- FINAL table rows: 1419

Final table

Show 5 entries Search:

row ID	row m/z	row retention time	row ion mobility	row ion mobility unit	row CCS	correlation group ID	annotation network number
1	108	148.0277472	0.3691534			31	0
2	181	169.0462455	0.369998			31	0
4	231	216.9512638	0.37140			31	0
6	145	244.9475245	0.37269926			31	0
8	189	217.98956	0.375543			31	23

Showing 1 to 5 of 1,419 entries Previous 1 2 3 4 5 ... 284 Next

3 Matches in filtered .mgf

by Plyushchenko I.V. | 2026 | [GitHub](#) | [Cite](#)



Tab “About” Overview

The screenshot displays the MetaboCensoR web application interface. At the top, there is a navigation bar with several tabs: 'Upload Data', 'Blank Filters', 'MS Filters', 'QC Filters', 'Peak Filters', 'Final Summary', and 'About'. The 'About' tab is currently selected. Below the navigation bar, the main content area is divided into two sections. On the left, there is a 'How to:' section with a dropdown menu titled 'Select a section:'. The dropdown menu lists the following options: 'Quick Start', 'Upload Data', 'Blank Filters', 'MS Filters', 'QC Filters', 'Peak Filters', 'Final Summary', 'Troubleshooting', and 'Project Details'. A red arrow points from the 'Troubleshooting' option in the dropdown menu to a yellow text box on the right. The yellow text box contains the text: 'Navigate through the sections to get more information'. On the right side of the main content area, there is a 'Troubleshooting' section with a list of bullet points:

- **No sample columns:** your sample keywords don't match sample column names.
- **Parsing error:** check selected Data table type and compatibility with examples of accessible Peak Table data formats in the [GitHub](#) repository.
- **Disconnected from the server:** could be caused by memory limit, try to compile and run it locally. Instructions in the [GitHub](#) repository.

At the bottom of the page, there is a footer with the text: 'by Plyushchenko I.V. | 2026 | [GitHub](#) | [Cite](#)'.