



MetaboCensoR Tutorial

2026-06-06

Tutorial Edition 1.2

Software Release 1.2.x



Contents

General Information.....	3
Running & Installation.....	4
Input Requirements.....	4
Overview of Interface	5
Implementation Details.....	6
Quick Start	6
Tab “Upload Data”: Control Panel.....	7
Tab “Upload Data”: Display Area.....	8
Tab “Blank Filters”: Control Panel	9
Tab “Blank Filters”: Display Area	10
Tab “MS Filters”: Control Panel.....	11
Tab “MS Filters”: Display Area	12
Tab “QC Filters”: Control Panel	13
Tab “QC Filters”: Display Area	14
Tab “Peak Filters”: Control Panel.....	15
Tab “Peak Filters”: Display Area	17
Tab “Final Summary”: Control Panel	15
Tab “Final Summary”: Display Area	18
Tab “About” Overview	19



General Information

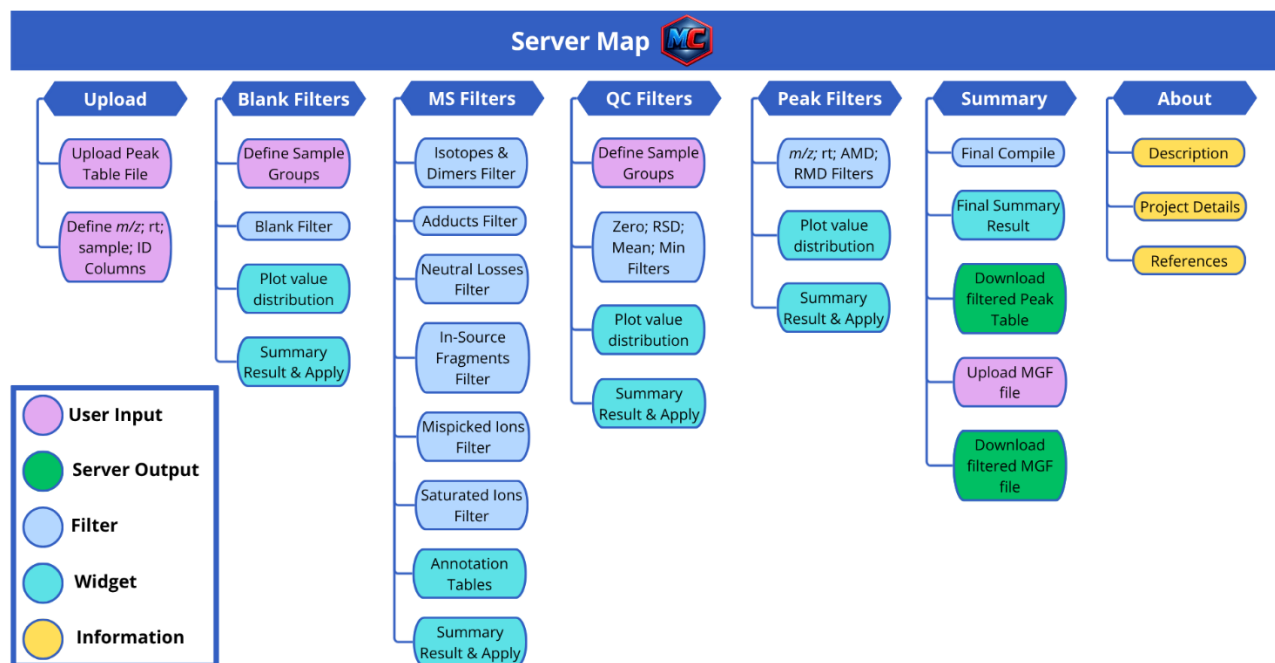
This tutorial aims at guiding the user to process data LC-MS data through Shiny GUI interface of MetaboCensoR. MetaboCensoR is an input-versatile Shiny application and local R package for analyte-centric peak table filtering. The workflow integrates four complementary modules for blank filtering, redundant ion-species filtering, quality-control filtering, and peak-based filtering. MetaboCensoR also provides interactive threshold optimization, exportable annotation tables, and synchronized filtering of associated .mgf files.

An overview of MetaboCensoR is presented below. The workflow is organized into six functional tabs that collectively support peak table curation, parameter optimization, and export of filtered results. The output of the app is a filtered peak table and, optionally, an MGF spectra file. Only the peak table upload is required for full workflow execution and final output generation. All other tabs are optional.

License: GPL-3.0

Contact: plyushchenko.ivan@gmail.com

Project Repository: github.com/plyush1993/MetaboCensoR





Running & Installation

Shiny Web Deployment is available by the link:

plyush1993.shinyapps.io/metabocensor/

Local installation of R package:

```
if (!require("BiocManager", quietly = TRUE)) {  
  install.packages("BiocManager")  
}  
if (!require("remotes", quietly = TRUE)) {  
  install.packages("remotes")  
}  
remotes::install_github("plyush1993/MetaboCensor", INSTALL_opts = "--no-multiarch")
```

or

```
if (!requireNamespace("pak", quietly = TRUE)) install.packages("pak")  
pak::pak("plyush1993/MetaboCensor")
```

Run Locally:

```
MetaboCensor::run_metabocensor()
```

Input Requirements

The application accepts peak tables in **.csv** format and MS1/MS2 mass spectral data in **.mgf** format, both of which can be directly exported from dedicated metabolomics processing software (xcms, MZmine, MS-DIAL). Make sure that your input peak table is aligned with data formats. Examples of acceptable formats are available here: [Input Examples](#).

Important note

Although the web version has been successfully tested on datasets with over 10,000 peaks, the Shiny server has limited memory capacity. Computationally intensive steps (like MS filters) or large MGF files (> 3 MB) may cause the app to stack or disconnect. For heavy processing, you should run the application locally.



Overview of Interface

The screenshot shows the user interface with three main panels highlighted by red boxes and numbers:

- 1 (Sidebar Panel):** Contains the 'Upload' section with a 'Data table type' dropdown (set to 'mzMine'), an 'Upload Peak Table (*.csv)' section with a 'Browse...' button, and 'Global parsing settings' including 'Feature ID (auto-generated)', 'm/z column', 'rt column', and 'How to define sample columns?' options.
- 2 (Main Panel):** Displays summary metrics: 'TOTAL SAMPLES DETECTED' (9) and 'TOTAL FEATURES (PEAKS)' (6,167). Below this is a 'Dataset: Example' section with a 'Show 6 entries' dropdown and a search bar. A table of 6 entries is shown with columns: row ID, row m/z, row retention time, row ion mobility, row ion mobility unit, row CCS, correlation group ID, annotation network number, best ion, auto MS2 verify, identified by n, partners, and neutral M mass.
- 3 (Navigation Bar):** A horizontal tab menu at the top with buttons: 'Upload Data', 'Blank Filters', 'MS Filters', 'QC Filters', 'Peak Filters', 'Final Summary', and 'About'.

GUI is organized into 3 main panels:

1. The Sidebar Panel (Control Panel)

- This is the primary area for user interaction. It houses the input controls where users upload data, configure global settings, define parameters, and execute filtering steps.

2. The Main Panel (Display Area)

- This is where the application renders the output and displays interactive data tables, statistical plots, and provides summary metrics.

3. The Navigation Bar (Tab Menu)

- This allows users to navigate linearly through the sequential steps of the data processing pipeline, from data upload to final summary and export.

Contextual assistance is embedded throughout the interface. Users can access detailed parameter descriptions and recommendations by interacting with the tooltip icons.

On the following pages, yellow numbered circles  indicate the recommended operational order within each interface tab in *Control Panel*.

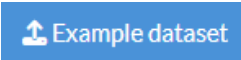
Red numbered squares  describe widgets in *Display Area*.



Implementation Details

All filter groups are applied sequentially in the order shown in the tab panel: Blank Filters → MS Filters → QC Filters → Peak Filters. Changes made on later tabs do not affect the results of previous tabs. To change the order of applied filters (e.g., 'Peak Filters' before 'Blank Filters'), run the application sequentially several times, using the output peak table from one run as the input for the next.

Quick Start

Press button  in “Upload data” tab to upload attached example dataset and get access to the MetaboCensoR functionality. If you use this example, run the Blank filter first to save memory for further steps. You can keep all parameters by default.

Briefly, it is the LC-MS plant profiling dataset, which is described in [Houriet et al.](#) Briefly, methanol extracts from the plant ashwagandha [*Withania somnifera* (L.) Dunal] together with blanks were analyzed on RP column in DDA positive mode on Thermo Q-Exactive Plus Orbitrap. Raw mzML files were then processed in mzMine in default pre-settings for UPLC-DDA to obtain resulted peak table.



Tab “Upload Data”: Control Panel

The screenshot shows the 'Upload Data' tab of the MetaboCensoR application. The interface includes a top navigation bar with tabs: Upload Data, Blank Filters, MS Filters, QC Filters, Peak Filters, Final Summary, and About. The main content area is divided into two columns. The left column contains the 'Upload' section with a 'Data table type' dropdown set to 'mzMine', a 'Browse...' button for uploading a peak table, and 'Global parsing settings' for Feature ID, m/z, and retention time columns. The right column displays 'TOTAL SAMPLES DETECTED' (9) and 'TOTAL FEATURES (PEAKS)' (6,167), followed by a 'Dataset: Example' section with a table of 6 entries. A yellow callout box on the right lists five steps: 1. Pick the right table type and upload; 2. Confirm the three core columns; 3. Specify sample intensity columns; 4. Clear dataset; 5. Example dataset shortcut. A yellow bracket on the left side of the 'Upload' section groups the 'Data table type' dropdown, the 'Browse...' button, and the 'Global parsing settings' section.

Upload

Data table type: ?
mzMine

Upload Peak Table (*.csv)
Browse... No file selected

Global parsing settings

Feature ID (auto-generated):
feature_id

m/z column:
row m/z

rt column:
row retention time

How to define sample columns?
☐ Auto-detect numeric sample columns
☒ By keyword match
☐ Pick columns manually

Sample column keywords (pick/add multiple):
.mzML .mzXML

Clear dataset

Example dataset

TOTAL SAMPLES DETECTED
9

TOTAL FEATURES (PEAKS)
6,167

Dataset: Example

Show 6 entries

Search:

row ID	row m/z	row retention time	row ion mobility	row ion mobility	row CCS	correlation group ID	annotation network number
1							
2							
3							
4							
5							
6							

Showing

next

- 1. Pick the right table type and upload**
MetaboCensoR supports mzMine, xcms, MS-DIAL, and Default input tables. If the type is wrong, the app throws a parsing error.
- 2. Confirm the three core columns**
Use the generated selectors to check Feature ID (by default auto-generates numeric vector), m/z, and rt. We recommend to start with default column names.
- 3. Specify sample intensity columns**
Use auto-detect, keyword matching, or manual column selection. Keyword matching is the most transparent because you can inspect the exact suffixes/prefixes used in your export.
- 4. Clear dataset**
Reset all calculations.
- 5. Example dataset shortcut**
The Example dataset button loads the built-in example.

by Plyushchenko I.V. | 2026 | GitHub | Cite



Tab “Upload Data”: Display Area



[Upload Data](#)[Blank Filters](#)[MS Filters](#)[QC Filters](#)[Peak Filters](#)[Final Summary](#)[About](#)

Upload

Data table type: ?
mzMine

Upload Peak Table (*.csv)
[Browse...](#) No file selected

Global parsing settings

Feature ID (auto-generated):
feature_id

m/z column:
row m/z

rt column:
row retention time

How to define sample columns?
☐ Auto-detect numeric sample columns
☒ By keyword match
☐ Pick columns manually

Sample column keywords (pick/add multiple):
.mzML .mzXML

[Clear dataset](#)

[Example dataset](#) ?

TOTAL SAMPLES DETECTED
9

1 Statistics of uploaded dataset

Dataset: Example
Show 6 entries Search:

row ID	row m/z	row retention time	row ion mobility	row ion mobility unit	row CCS	correlation group ID	annotation network number
1	108	148.0277472	0.3691534			31	0
2	181	169.0462455	0.3699988			31	
3	175	128.0197664	0.37076914			31	20
4	231	216.951263					
5	43	234.961351	0.3725722			31	14
6	145	244.9475245	0.37269926			31	

Showing 1 to 6 of 6,167 entries [Previous](#) [1](#) [2](#) [3](#) [4](#) [5](#) ... [1,028](#) [Next](#)

by Plyushchenko I.V. | 2026 | [GitHub](#) | [Cite](#)



Tab “Blank Filters”: Control Panel

The screenshot shows the 'Blank Filters' tab in the MetaboAnalyst interface. The top navigation bar includes 'Upload Data', 'Blank Filters' (selected), 'MS Filters', 'QC Filters', 'Peak Filters', 'Final Summary', and 'About'. The left sidebar contains the 'Sample Label' section with three options: 'From sample names (by token)', 'From custom CSV (one column)', and 'Manual editable table'. The 'Manual editable table' option is selected. Below this is a 'Blank / Media Removal' section with an 'Enable Blank filter' toggle, a 'Choose blank/media group(s)' dropdown, and a 'Filtering mode' section with two options: 'Blank < Cutoff * Sample Max' (selected) and 'Drop any feature detected in Blank'. A 'Cutoff (fraction of max experimental mean)' input field is set to 0.1. There are 'Apply' and 'Clear' buttons. The main area shows an 'Input table' with 6 rows and 6 columns. Below the input table is a 'Blank to' histogram. The 'Output table' shows 6 rows and 6 columns. At the bottom, a 'Filtering results — Step 1 (Blank Filters)' summary table is displayed.

row	ID	1	2	3	4
1	108	148.0277472	0.3691534	31	0
2	181	169.0462455	0.3699988	31	
4	231	216.9512638	0.371403	31	
5	43	234.961351	0.3725722	31	14
6	145	244.9475245	0.37269926	31	

Metric	Value
Features before	6167
Removed by Blank Filters	3570
Features after	2597

1. Sample Labels — how to specify groups

Choose a label source: from sample names or custom CSV (one column, no header).

Define token and index for extracting labels from sample names.

Example — parse from sample names

Sample name: Sample_WT_1_R-C1_1_5037.mzML Peak area

Separator: _

Token index: 2

Resulting label: WT

From custom CSV:

Upload a one-column CSV without header. Label order must match sample column order.

Manual editable table:

The editable label table is filled from the current token-based labels. Double-click cells in the Label column to manually correct group names.

2. Show resulted Labels in Display Area

3. Set Up Blank Filters

Enable filter.

Specify blank groups.

Choose between filtering modes:

- Blank < Cutoff * Sample Max
- Drop any feature detected in Blank

4. Plot distribution of filtered peaks

5. Apply or Discard Filtering Step



Tab “Blank Filters”: Display Area

Upload Data

Blank Filters

MS Filters

QC Filters

Peak Filters

Final Summary

About

Sample Label

Label source:
☐ From sample names (by token)
☐ From custom CSV (one column, no header)
☒ Manual editable table

Fill editable table from current token labels

☒ Show labels table

Blank / Media Removal

Enable Blank filter ☒

Choose blank/media group(s):
Pick one or more (e.g., Media, Blank)

Filtering mode:
☒ Blank < Cutoff * Sample Max
☐ Drop any feature detected in Blank

Cutoff (fraction of max experimental mean):
0.1

Plot values distribution

☒ Show Blank plot panel

Apply

Clear

Input table

Show 5 entries

Search:

row ID	mz	rt	row ion mobility	row ion mobility unit	row CCS	correlation group ID	annotation network number
1	108	148.0277472	0.3691534			31	0
2	181	169.0462455	0.3699988			31	
3	175	128.0197664	0.37076914			31	20
4	231	216.9512638				31	
5	43	234.961351				31	14

Showing 1 to 5 of 6,167 entries

Previous12345...1,234Next

Check matching Sample Names and Group Labels:

Show 6 entries

Sample	Label
Orbi_Blank_A.mzML Peak area	Blank
Orbi_Blank_B.mzML Peak area	Blank
Orbi_Blank_C.mzML Peak area	y
Orbi_Sample_A.mzML Peak area	Sample
Orbi_Sample_B.mzML Peak area	Sample
Orbi_Sample_C.mzML Peak area	Sample

Showing 1 to 6 of 9 entries

Previous12Next

Check Blank Ratio Distribution:

Blank to Sample Ratio Distribution

3
Distribution plot

Output table — after Blank Filters

Show 5 entries

Search:

row ID	mz	rt	row ion mobility	row ion mobility unit	row CCS	correlation group ID	annotation network number
1	108	148.0277472	0.3691534			31	0
2	181	169.0462455	0.3699988			31	
4	231	216.9512638	0.371403			31	
5	43	234.961351	0.372				14
6	145	244.9475245	0.372				

Showing 1 to 5 of 2,597 entries

Previous12345...520Next

Filtering results — Step 1 (Blank Filters)

Metric	Value
Features before	6167
Removed by Blank Filters	3570
Features after	2597

5
Summary



Tab “MS Filters”: Control Panel

A) Isotopes & Dimers

Enable isotope/dimer collapse ☒

Correlation Z: 0.8

rt tolerance (min): 0.01

m/z tolerance type: ☒ Da ☐ ppm

m/z tolerance (Da): 0.005

Max charge (z_max): 3

Max isotope order (n): 1

Dimer-series max n_d: 3

B) Adduct Families

Enable adduct-family collapse ☒

Polarity: ☒ positive ☐ negative

Adducts CSV (optional):

Correlation Z: 0.8

rt tolerance (min): 0.005

m/z tolerance type: ☒ Da ☐ ppm

m/z tolerance (Da): 0.005

Minimum neutral mass: 50

☒ Strict RT split inside clusters

C) Neutral Losses

Enable neutral loss collapse ☐

D) In-Source Fragments

Enable in-source fragment collapse ☐

E) Mispicked Ions

Enable mispicked ions collapse ☐

F) Saturated Ions

Enable ringing artifact collapse ☐

☒ Show isotope/dimer groups

☐ Show adduct families

☐ Show adduct statistics

☐ Show neutral loss hits

☐ Show fragment clusters

☐ Show mispicked ions groups

☐ Show ringing artifact families

Input table

Show 5 entries

Search:

row ID	mz	rt	row ion mobility	row ion unit	row CCS	correlation group ID	annotation number
1	108	148.0277472	0.3691534			31	0
2	181	169.0462455	0.3699988			31	
3	175	128.0197664	0.37076914			31	20
4	151	125.12638	0.371403			31	
5	1351	0.37257					

Output table - after MS

Show 5 entries

row ID	mz	correlation group ID	
1	108	148.0277472	0.3691534
3	175	128.0197664	0.37076914
7	177	151.0356683	0.37076914
10	171	110.0093424	0.3785
11	141	223.989142	0.379040

Showing 1 to 5 of 4,258 entries

Isotopes table:

Feature	mz	correlation group ID
981.5019	5.5	
982.5019	5.5	
977.9383	9.2	
977.4369	9.2	
972.4807	9.2	
972.9820	9.2	
895.9259	4.4	
896.4304	4.4	
883.8535	4.4	
884.3593	4.4	

Showing 1 to 5 of 6,167 entries

1. Enable & Set Up Isotopes Filtering

Need parameter help? Click the tooltips or read the [About](#) tab for details. **Pro-Tip:** Start with the default settings.

2. Enable & Set Up Adducts Filtering

Need parameter help? Click the tooltips or read the [About](#) tab for details. **Pro-Tip:** Start with the default settings and check your polarity.

3. Enable & Set Up Neutral Losses Filtering

Need parameter help? Click the tooltips or read the [About](#) tab for details. **Pro-Tip:** Start with the default settings and check your polarity.

4. Enable & Set Up In-Source Fragments Filtering

Need parameter help? Click the tooltips or read the [About](#) tab for details. **Pro-Tip:** Start with the default settings.

4. Enable & Set Up Mispicked Ions Filtering

Need parameter help? Click the tooltips or read the [About](#) tab for details.

6. Enable & Set Up Saturated Ions Filtering

Need parameter help? Click the tooltips or read the [About](#) tab for details.

7. Print & Download annotation tables

8. Apply or Discard Filtering Step



12



Tab “QC Filters”: Control Panel

The screenshot shows the QC Filters control panel with the following components:

- 1. Sample Labels**: A section for specifying label sources, including "Use labels from Blank Tab", "From sample names (by token)", "From custom CSV (one column, no header)", and "Manual editable table".
- 2. QC value filters**: A section for enabling filters (Zeros, Mean, RSD, Min) and setting evaluation modes (Groupwise, Pooled).
- 3. Zeros metric**: A section for setting the Zeros metric (Count or Percent) and the Zeros threshold (count or %).
- 4. Mean/RSD/Min thresholds**: A section for setting the Mean z threshold, RSD z threshold, and Min z threshold.
- 5. Plot values distribution**: A section for plotting the distribution of selected values, with a checkbox for "Show QC plot panel".
- 6. Output table**: A table showing the results of the QC filters, including row ID, mz, rt, and row ion mobility.
- 7. Filtering results**: A table showing the results of the filtering process, including the number of features before and after each filter step.

Metric	Value
Features before	2597
Removed by Zeros	5
Removed by Mean	1453
Removed by RSD	110
Removed by Min	595
Total removed	1485
Features after	1112

1. Sample Labels — how to specify groups

Choose a label source: from *Blank Filters* tab, sample names, custom CSV (one column, no header), or by manual editing.

See details for extracting Labels on the Page **Tab “Blank Filters”: Control Panel**.

2. Show resulted Labels in Display Area

3. Set Up QC Values Filters

Select desired values for filtering:

Zeros, Mean, RSD, or Min.

Choose between filtering modes:

- Groupwise (ANY selected group passes)
- Groupwise (EVERY selected group passes)
- POOLED across selected groups

Specify groups used for evaluation.

4. Set Up Zero Values Filters

Specify Zero metric (counts or %) and threshold.

5. Set Up Mean/RSD/Min Values Filters

Specify thresholds for corresponding value.

6. Plot distribution of selected values

7. Apply or Discard Filtering Step



Tab “QC Filters”: Display Area

MC

MetaboAnalyst 5.0

Upload Data

Blank Filters

MS Filters

QC Filters

Peak Filters

Final Summary

About

Sample Label

Label source:

Use labels from Blank Tab

From sample names (by token)

From custom CSV (one column, no header)

Manual editable table

Show labels table

QC value filters

Enable filters:

Zeros

Mean

RSD

Min

Evaluation mode:

Groupwise (ANY selected group passes)

Groupwise (EVERY selected group passes)

POOLED across selected groups

Group(s) to evaluate:

Empty - all groups

Zeros metric:

Count

Percent

Zeros $\leq$ threshold (count or %)

0

Mean $\geq$ threshold

5000

RSD $\leq$ threshold

30

Min $\geq$ threshold

1000

Plot values distribution

Show QC plot panel

Apply

Clear

Input table

Show 5 entries

Search:

row ID	mz	rt	row ion mobility unit	row CCS	correlation group ID	annotation number
1	108	148.0277472	0.3691534		31	0
2	181	169.0462455	0.3699988		31	
4	231	216.95			31	
5	43	234.9			31	14
6	145	244.9475245	0.37269926		31	

Showing 1 to 5 of 2,597 entries

Previous12345...520Next

Check matching Sample Names and Group Labels:

Show 6 entries

Sample	Label
Orbi_Blank_A.mzML Peak area	Blank
Orbi_Blank_B.mzML Peak area	Blank
Orbi_Blank_C.mzML Peak area	Blank
Orbi_Sample_A.mzML Peak area	Sample

Check Values Distribution:

ZEROS [Zero values (count)] - Groupwise (ANY) | Filtered: 5 of 2597 (0.2%)

MEAN [Mean intensity (log10)] - Groupwise (ANY) | Filtered: 1453 of 2597 (55.9%)

MIN [Minimum intensity (log10)] - Groupwise (ANY) | Filtered: 595 of 2597 (22.9%)

Distribution plots

Output table — after QC Filters

Show 5 entries

Search:

row ID	mz	rt	row ion mobility unit	row CCS	correlation group ID	annotation number
1	108	148.0277472	0.3691534		31	0
2	181	169.0462455	0.3699988		31	
4	231	216.9			31	
5	43	234			31	14
6	145	244.9475245	0.37269926		31	

Showing 1 to 5 of 1,112 entries

Previous12345...223Next

Filtering results — Step 3 (QC Filters)

Metric	Value
Features before	2597
Removed by Zeros	5
Removed by Mean	1453
Removed by RSD	110
Removed by Min	5
Total removed	35
Features after	1112

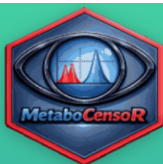
Note: a feature may be counted in multiple filter categories.

by ProtonShine MS (1.65x2) v. 1.3.0 | Publication

14



Tab “Final Summary”: Control Panel

[Upload Data](#)[Blank Filters](#)[MS Filters](#)[QC Filters](#)[Peak Filters](#)[Final Summary](#)[About](#)

Compile final output

[Compile output](#)[Final peak table \(.csv\)](#)

MGF Filtering

Filter MGF

Upload MGF file

[Browse...](#)

mgf_gnps_test.mgf

Upload complete

Select ID column in Table:

row ID

MGF Feature ID field:

acquisitionNum

[Filter MGF](#)[Clear MGF](#)

Success! Matched 2 spectra from your filtered table.

[Download Filtered MGF](#)

Filtering summary

- Dataset: Example
- Step 1 Blank Filters: APPLIED
- Step 2 MS Filters: APPLIED

1. Compile final filtered table

Clicking this button enables the table download and MGF filtering.

2. Download final table

3. Enable MGF filtering & upload .mgf file

4. Match peak ID between table and .mgf file

The feature identifier in the final peak table must correspond to the chosen spectrum identifier field inside the uploaded MGF. If no IDs match, the app reports: “No matching spectra were found in the MGF file”.

We recommend using the default columns unless you have specific requirements.

5. Filter or clear uploaded .mgf file

6. Download filtered .mgf file



Tab “Peak Filters”: Control Panel

Peak Filters

Enable filters: ☒ m/z ☒ rt ☐ RMD ☐ AMD

☒ m/z ≥ 200 ☐ m/z ≤

☒ rt ≥ (min) 1 ☐ rt ≤ (min)

☒ Plot values distribution

☒ Show Peak plot panel

Input table

Show 5 entries Search:

row ID	m/z	rt	row ion mobility	row ion mobility unit	row CCS	correlation group ID	annotation network number
1	108	148.0277472	0.3691534			31	0
2	169	0.462455	0.3699988			31	
3	216	9512638	0.371403			31	
4	234	961351	0.3725722			31	14
5	244	9475245	0.37269926			31	

Showing 1 to 5 of 1,068 entries Previous 1 2 3 4 5 ... 214 Next

Plot values distribution

Filtered m/z: 157
Filtered rt: 316
Total Filtered: 369

1. Set Up Peak Values Filters

Select desired values for filtering:
m/z, *rt*, or mass defect (RMD, or AMD).

2. Set Up *m/z*/*rt*/AMD/RMD Values Filters

Specify threshold ranges for corresponding value.

3. Plot distribution of selected values

4. Apply or Discard Filtering Step

Output table — after Peak Filters

Show 5 entries Search:

row ID	m/z	rt	row ion mobility	row ion mobility unit	row CCS	correlation group ID	annotation network number
818	2815	336.074296	1.0103563			58	489
819	2820	314.0921816	1.0123757			58	489
825	2879	291.0843679	1.0544833			60	495
829	2893	275.1104506	1.0585195			60	495
830	2862	253.1291676	1.0596385			60	495

Showing 1 to 5 of 699 entries Previous 1 2 3 4 5 ... 140 Next

Filtering results — Step 4 (Peak Filters)

Metric	Value
Features before	1068
Removed by <i>m/z</i>	157
Removed by <i>rt</i>	316
Removed by RMD	0
Removed by AMD	0
Total removed	369



Tab “Peak Filters”: Display Area

MC

MetaboAnalyst 5.0

Upload Data

Blank Filters

MS Filters

QC Filters

Peak Filters

Final Summary

About

Peak Filters

Enable filters: 1

m/z

rt

RMD

AMD

☒ m/z ≥

☐ m/z ≤

☒ rt ≥ (min)

☐ rt ≤ (min)

200

1

Plot values distribution

☒ Show Peak plot panel

Apply

Clear

Input table

row ID	mz	rt	row ion mobility	row CCS	correlation group ID	annotatic network number
1	108	148.0277472	0.3691534		31	0
2	181	169.0462455	0.3699988		31	
4	231	216.9512638	0.372			14
5	43	234.961351	0.372			
6	145	244.9475245	0.37269926		31	

Showing 1 to 5 of 1,068 entries

Distribution plot

m/z vs RT

Filtered m/z: 157
Filtered rt: 316
Total Filtered: 369

Output table

row ID	mz	rt	row ion mobility unit	row CCS	correlation group ID	annotatic network number
818	2815	336.07429				489
819	2820	314.092181				489
825	2879	291.0843679	1.0544833		60	495
829	2893	275.1104506	1.0585195		60	495
830	2862	253.1291676	1.0596385		60	495

Showing 1 to 5 of 699 entries

Summary

Metric	Value
Features before	1068
Removed by m/z	157
Removed by rt	316
Removed by RMD	0
Removed by AMD	0
Total removed	369
Features after	699

Note: a feature may be counted in multiple filter categories.



Tab “Final Summary”: Display Area



[Upload Data](#)[Blank Filters](#)[MS Filters](#)[QC Filters](#)[Peak Filters](#)[Final Summary](#)[About](#)

Compile final output

Compile output ?

Final peak table (.csv)

MGF Filtering

Filter MGF ☒ ?

Upload MGF file

Browse...

mgf_gnps_toy.mgf

Upload complete

Select ID column in Table:

row ID

MGF Feature ID field:

acquisitionNum

Filter MGF

Clear MGF

Success! Matched 2 spectra from your filtered table.

Download Filtered MGF

Filtering summary

- Dataset: Example
- Step 1 Blank Filters: APPLIED
- Step 2 MS Filters: APPLIED
- Step 3 QC Filters: not applied
- Step 4 Peak Filters: not applied
- FINAL table rows: 1419

Final table

Show 5 entries Search:

row ID	row m/z	row retention time	row ion mobility	row ion mobility unit	row CCS	correlation group ID	annotation network number
1	108	148.0277472	0.3691534			31	0
2	181	169.0462455	0.369998			31	0
4	231	216.9512638	0.37140			31	0
6	145	244.9475245	0.37269926			31	0
8	189	217.98956	0.375543			31	23

Showing 1 to 5 of 1,419 entries Previous 1 2 3 4 5 ... 284 Next

3 Matches in filtered .mgf

by Plyushchenko I.V. | 2026 | [GitHub](#) | [Cite](#)



Tab “About” Overview

The screenshot displays the MetaboCensoR web application interface. At the top, there is a navigation bar with the following tabs: Upload Data, Blank Filters, MS Filters, QC Filters, Peak Filters, Final Summary, and About. The 'About' tab is currently selected. On the left side, there is a 'How to:' dropdown menu with the following options: Quick Start, Upload Data, Blank Filters, MS Filters, QC Filters, Peak Filters, Final Summary, Troubleshooting, and Project Details. A red arrow points from the 'Troubleshooting' option in the dropdown menu to a yellow text box that reads: 'Navigate through the sections to get more information'. The main content area of the 'About' tab is titled 'Troubleshooting' and contains three bullet points:

- **No sample columns:** your sample keywords don't match sample column names.
- **Parsing error:** check selected Data table type and compatibility with examples of accessible Peak Table data formats in the [GitHub](#) repository.
- **Disconnected from the server:** could be caused by memory limit, try to compile and run it locally. Instructions in the [GitHub](#) repository.

At the bottom of the page, there is a footer that reads: 'by Plyushchenko I.V. | 2026 | [GitHub](#) | [Cite](#)'.